

PRECISE ANATOMICAL DISSECTION

Maximize your surgical technique and optimize outcomes

C2Dx® HEMOSTATIX®
Thermal Scalpel



THERMAL SCALPEL SYSTEM

THE ONLY SURGICAL
SCALPEL THAT SEALS
VESSELS AS IT CUTS

Experience a clean surgical field with exceptional visualization for the safe dissection of nerves and critical anatomy. No muscle excitation, no electric current.

A NEW PARADIGM FOR

FINE DISSECTION
NERVE PRESERVATION
TISSUE INTEGRITY



SURGICALLY SHARP

Precise excision of metastatic disease, unparalleled neural protection and preservation of vital anatomical structures.



THERMAL HEMOSTASIS

Heat from the scalpel blade is transferred to incised tissue within a narrow margin and no undesirable spread.



MAXIMUM CONTROL

Control the rate, depth, and angle of dissection along with blade temperature - while receiving accurate visual and tactile feedback.



SMOKE -FREE

Minimize capillary oozing and blood loss without electroconductivity. Cut and cauterize without creating harmful surgical smoke.

*To know it and to use it, is to
love it and then need it.*

— Mark Ochs, DMD, MD

See for yourself




40 Years of Clinical History Meets Leading Edge Technology

Printed circuit embedded
onto the blade surface

Blade hub with
magnetic latch


Energize/
De-energize
switch

COAG (max blade temp) button

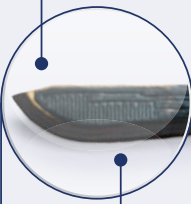

Electro-chemically honed blade edge
and state-of-the-art non-stick coating

Intuitive temperature controls

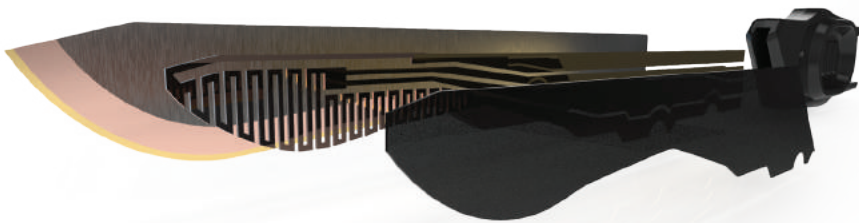

THE SCALPEL BLADE

A 5-LAYER LAMINATE COMPOSITION
OF COPPER AND STAINLESS STEEL

Advanced micro-electric circuitry applied to the blade surface enables immediate heating to a maximum of 300 degrees Celsius. The patient is insulated from the electric current, no grounding pad is required, and no muscle excitation occurs.

PRODUCT ORDERING INFORMATION

70138400	Precision 8400 Controller
70139050	Scalpel Handle for Precision 8400 Controller (1)
70239050	Scalpel Handle for Precision 8400 Controller (6)
70235810	#10 Scalpel Blade for Precision 8400 (10 per box)
70135810	#10 Scalpel Blade for Precision 8400 (24 per box)
70235812	#12 Scalpel Blade for Precision 8400 (10 per box)
70235815	#15 Scalpel Blade for Precision 8400 (10 per box)
70135815	#15 Scalpel Blade for Precision 8400 (24 per box)

