

Treat obstructive salivary gland disorders.

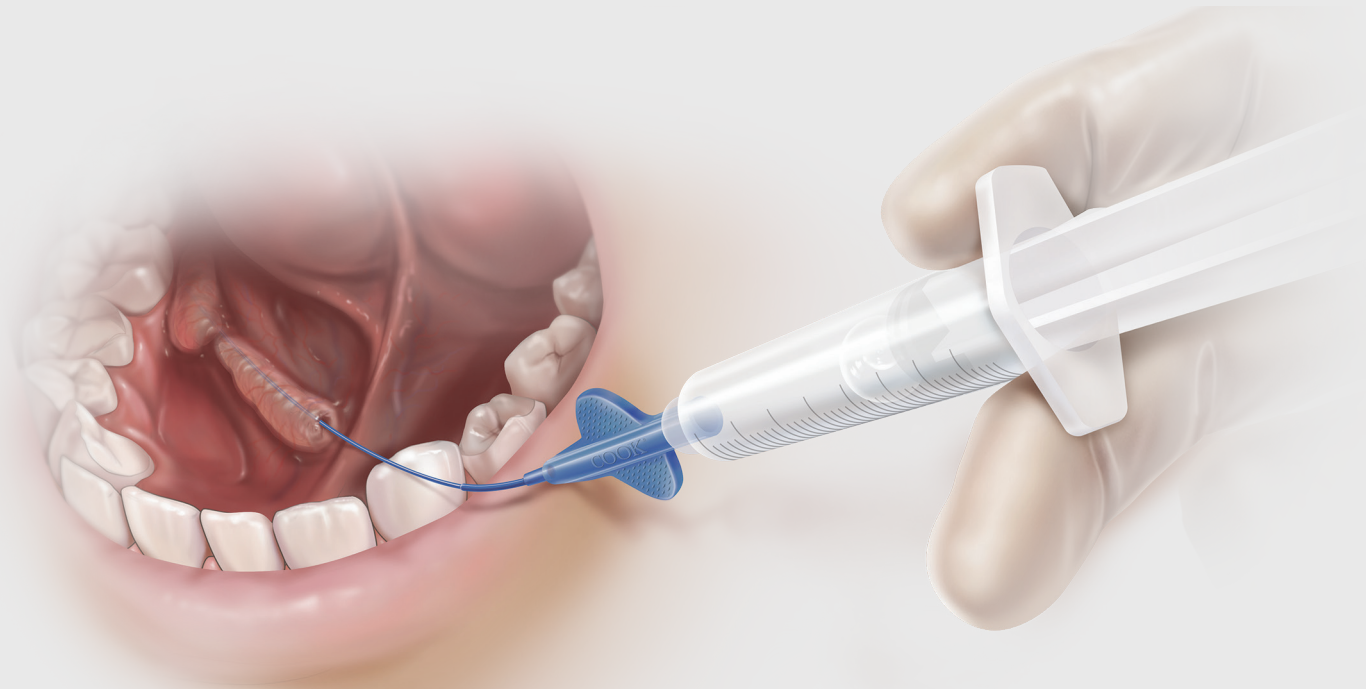


Illustration by Lisa Clark

SialoCath™

SALIVARY DUCT CATHETER

Scope-free irrigation


The catheter allows for irrigation without the use of a scope as a first-line, in-office treatment approach.

Atraumatic design

The soft tip, flexible catheter provides a gentle treatment option for both in-office and surgical procedures.

Salivary use

The catheter is specifically indicated for salivary duct irrigation.

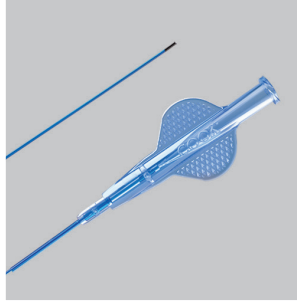

SialoCath™


SALIVARY DUCT CATHETER

Used to aid in the removal and treatment of salivary duct stones and other duct conditions through flushing and/or fluid irrigation of the salivary ducts.

This flexible, single-lumen catheter can be inserted into the salivary duct either over a wire guide or through an appropriately sized introducer sheath.

Order Number	Reference Part Number	Catheter OD mm/inch	Catheter ID mm/inch	Compatible Wire Guide Diameter inch	Compatible Catheter Length cm	Tip Configuration
G35166	SDC-2.2-20	.74/.029	.53/.021	.018	15	Straight


Related products

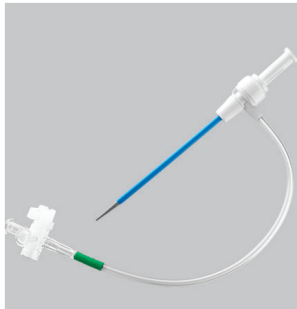

Kolenda Salivary Access Introducer Set

Used for the introduction of diagnostic or interventional devices in the submandibular salivary duct.

Order Number	Reference Part Number	Introducer ID mm/inch	Wire Guide Diameter inch	Sheath Length cm
G24129	SDAS-6.0-18	2.2/.087	.018	5.5
G24132	SDAS-5.0-18	1.9/.074	.018	5.5

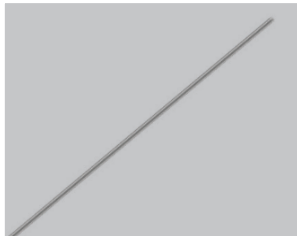

Wire guides

Used to gain and maintain access to the salivary duct.

Order Number	Reference Part Number	Wire Guide Diameter inch	Wire Guide Length inch	Tip Configuration
G02274	C-SF-15-50	.015	50	Straight

Some products or part numbers may not be available in all markets. Contact your local C2dx sales representative or Customer Service for details.



* Copyright© 2024 C2Dx
 * Printed in U.S.A.
 * 555 East Eliza St, Ste. A
 Schoolcraft, MI 49087 USA
 * t: 888 902 2239
 www.c2dx.com